

Rabbit MaxPab[®]

- A New Breed of Polyclonal Antibodies

Rabbit MaxPab[®] is a new collection of rabbit polyclonal antibodies generated from *full-length* proteins instead of peptide antigens for immunization. Rabbit MaxPab[®] antibody has the ability to recognize both linear and conformational protein epitopes for maximum antibody performance. In addition, each rabbit MaxPab[®] antibody is checked against a mammalian transfected lysate for demonstration of high antibody reactivity, sensitivity and specificity.

Rabbit MaxPab[®] Advantages

- Made from *full-length* proteins
- Multiple epitope recognition
- QC using mammalian transfected lysate
- Broad scope applications (WB, IF, IHC, IP, FC)
- Ideal capture antibody for ELISA pair study
- Higher antibody titer and affinity
- Larger sera batch for consistency and quality



Power Precision Performance


Transform your research with MaxPab[®]

Western Blot (Transfected Lysate)

Product Name : CCNE1 purified MaxPab[®] rabbit polyclonal antibody
Catalog # : H00000898-D01P

Western Blot analysis of CCNE1 expression in transfected 293T cell line by CCNE1 MaxPab[®] rabbit polyclonal antibody.


Lane 1: CCNE1 transfected lysate(47.1 KDa).
 Lane 2: Non-transfected lysate.

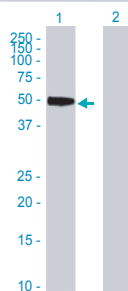

Western Blot (Tissue Lysate)

Product Name : AKR1A1 purified MaxPab[®] rabbit polyclonal antibody
Catalog # : H00010327-D01P


AKR1A1 MaxPab[®] rabbit polyclonal antibody. Western Blot analysis of AKR1A1 expression in human colon.

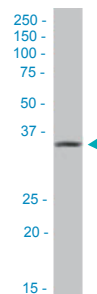

ELISA Antibody Pair

Product Name : Human FKBP4 Matched Antibody Pair
Catalog # : H00002288-AP31


- Capture antibody: affinity purified rabbit MaxPab[®] polyclonal anti-FKBP4
- Detection antibody: biotinylated affinity purified MaxPab[®] polyclonal anti-FKBP4

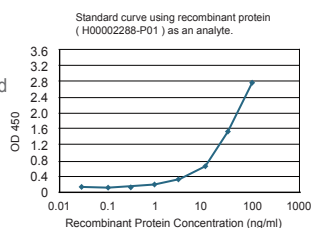

Immunohistochemistry

Product Name : GOT1 purified MaxPab[®] rabbit polyclonal antibody
Catalog # : H00002805-D01P

Immunoperoxidase of purified MaxPab[®] antibody to GOT1 on formalin-fixed paraffin-embedded human colon. [antibody concentration 3ug/ml]

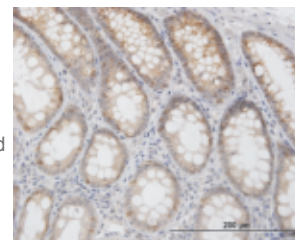

Immunofluorescence

Product Name : CALB1 purified MaxPab[®] rabbit polyclonal antibody
Catalog # : H00000793-D01P

Immunofluorescence of purified MaxPab[®] antibody to CALB1 (H00000793-D01P) on HeLa cell. [antibody concentration 10ug/ml]



Flow Cytometry

Product Name : FAS purified MaxPab[®] rabbit polyclonal antibody
Catalog # : H00000355-D01P

Surface staining of human peripheral blood cells with anti-human CD95 (FAS) FITC.

